
vetscan Imagyst®

Digital Cytology and Add-on Expert Review

USER GUIDE

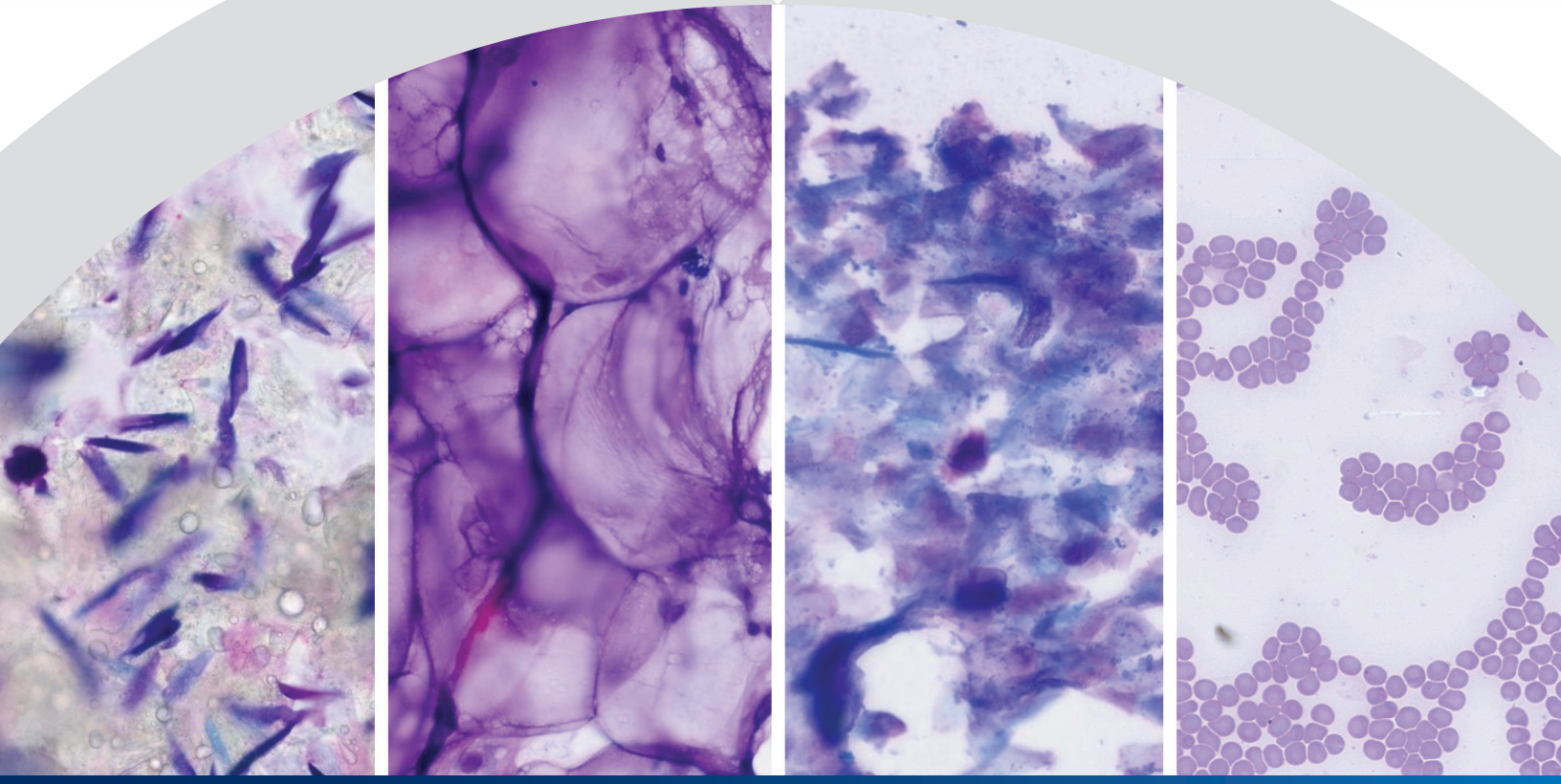

Getting started with Vetscan Imagyst for efficient Digital Cytology results and Add-on Expert review for Vetscan Imagyst AI tests*

vetscanimagyst.com

Efficient, expert cytology review


Vetscan Imagyst accelerates the process for expert clinical pathology review—simply prepare and scan slides and then submit for expert results

01 PREPARE CYTOLOGY SLIDE




- Prepare a cytology slide using industry best practices
- Place a drop of immersion oil on the sample and then place a 24 x 60 mm coverslip on the slide
- Include a label or handwritten note on the slide
- Keep the slides together and well organized

02 ADD NEW TEST

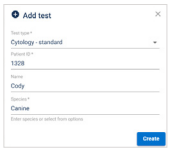

- Log in to Vetscan Imagyst
- Choose the correct sample and select **START**. Sample information will come prepopulated from any Fuse- or Hub-connected software system
- If no integration is available, select **Add Test (+)** and choose Standard or Urgent service
- Enter all required information and then select **CREATE**

03 COMPLETE PATIENT HISTORY

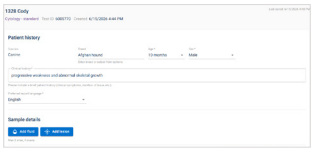


- Enter the patient's breed, age and gender
- The patient's species is already set
- Include a brief clinical history

04 ADD SAMPLE SITES




- Select **ADD FLUID** or **ADD LESION**
- Enter the required information
- Select **CREATE**
- You may create up to 2 different sample sites
- Select a sample card to open and edit a sample site
- To remove a sample site, select **DELETE (🗑)**

05 ADD SCANS

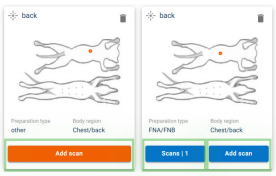


- Select **ADD A SCAN** or **SCANS** to open the scan window
- Select an available scanner
- Select **OPEN TRAY** and place cytology slide on the tray, locking it in place

06 CREATE SCAN AREA

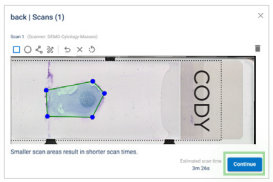

- Select **PREVIEW** and review the suggested scan area
- Accept the suggested scan area or create a custom scan area to reduce the scan time. Select **CONTINUE** when you are satisfied with the scan area
- You may submit up to 2 sample sites, with a maximum of 4 slides in total
- Select **DELETE (X)** to remove a scan

07 SUBMIT TEST

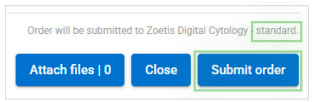


- To assist the pathologist with their diagnosis, attach any history documents (.pdf) or photos (.jpeg)
- Verify the cytology service displayed above **SUBMIT ORDER** button is correct
- Select **SUBMIT ORDER**

08 REVIEW TESTS



- When a test is ready to be reviewed, select **READY FOR REVIEW**
- Find the desired test and select **REVIEW**
- Review the findings
- The test will move to **TEST HISTORY**

INCOMPLETE TESTS

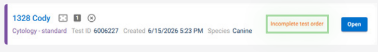

- An incomplete or unsubmitted test will display with an **Incomplete Test Order** status in the Inbox
- Select **OPEN** to complete all required fields
- Review steps 3-7

Contact your Zoetis representative for assistance.

zoetis

Efficient, expert cytology review


The Add-on Expert Review in Vetscan Imagyst allows you to quickly route AI-analyzed cases to a board-certified clinical pathologist for a second opinion.

01 PREPARE
SLIDE (Optional)




- For AI Urine Sediment - prepare a stained urine sediment smear
- For AI Masses – prepare stained slides from same site
- Refer to AI Urine Sediment Analysis or AI Masses Quick Start Guide for slide preparation best practices

02 ADD
NEW TEST

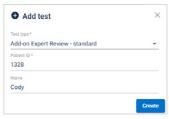

- Log in to Vetscan Imagyst
- Choose the correct sample and select **START**. Choose the correct AI test for which an expert review is required and click **SELECT TEST**. Sample information will come prepopulated from any Fuse or Hub connected software system
- If no integration is available, select **Add Test (+)** and choose Standard or Urgent service
- Enter all required information and then select **CREATE**

03 COMPLETE
PATIENT HISTORY

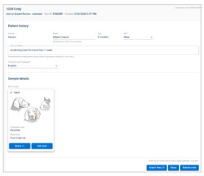


- Enter the patient's breed, age and gender
- The patient's species is already set
- Include a brief clinical history

04 PROVIDE
ADDITIONAL SLIDES (Optional)



- For AI Urine Sediment only 1 additional slide can be included
- For AI Masses up to 3 additional slides from the same site can be included for expert review
- An expert review for AI Blood Smear or AI Dermatology test does not require additional slides.

05 ADD
SCANS

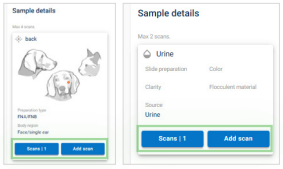

- Select **SCANS** or **ADD SCANS** to open the scan window
- Select an available scanner
- Select **OPEN TRAY** and place slide on the tray, locking it in place

06 CREATE
SCAN AREA



- Select **PREVIEW** and review the suggested scan area
- Accept the suggested scan area or create a custom scan area to reduce the scan time. Select **CONTINUE** when you are satisfied with the scan area
- Select **DELETE (X)** to remove a scan

07 SUBMIT
TEST




- To assist the pathologist with their diagnosis, attach any history documents (.pdf) or photos (.jpeg)
- Verify the expert review service displayed above **SUBMIT ORDER** button is correct
- Select **SUBMIT ORDER**

08 REVIEW
TESTS



- When a test is ready to be reviewed, select **READY FOR REVIEW**
- Find the desired test and select **REVIEW**
- Review the findings
- The test will move to **TEST HISTORY**

INCOMPLETE TESTS

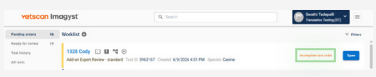

- An incomplete or unsubmitted test will display with an **Incomplete Test Order** status in the Inbox
- Select **OPEN** to complete all required fields
- Review steps 3-7

Contact your Zoetis representative for assistance.

*AI Blood Smear, AI Dermatology, AI Urine Sediment and AI Masses

All trademarks are the property of Zoetis Services LLC or a related company or licensor unless otherwise noted.
© 2026 Zoetis Services LLC. All rights reserved. VTS-00477 R3